

MAGNETOM Vida

Coil Quality Assurance

Serial Number: 175692

Coil Quality Assurance

Demographical Data					
Customer		System		Service	
Customer Name		Product Name	MAGNETOM Vida	Service Center Name	Customer Care Center
Customer ID	400-567654	Material Number	11060815	Street	Wheeler Road
Hospital Name	Beth Israel Deaconess	Serial Number	175692	Street Number	200 suite 303
Street	Brookline Ave	Equipment Number		Zip Code	
Street Number	330	Modality	MR	City	Boston
Zip Code	02215	Handed Over Date	6/28/2018	District	
City	Boston	Software Version	syngo MR XA11	Country	US
District	NorthEast			Phone	18008887436
Country	US			Fax	
Customer Phone Number	16179759913			E-Mail	uptime-services.team@siemens-healthineers.com

Summary		
Function	Start Time	Result
Coil Check - Flex Large 4 (30503)	5/14/2020 3:57:13 PM	 OK

1 Coil Check - Flex Large 4 (30503)

1.1 Summary


Summary	
	Status
ExamDb:SIEMENS/System/CoilQA/FlexLarge_4/099_FlexLarge_4_Tra	OK
Program Result	Success

1.2 Details

1.2.1 Protocol: ExamDb:SIEMENS/System/CoilQA/FlexLarge_4/099_FlexLarge_4_Tra

Coil Info		
Coil ID	Serial Number	Part Number
FlexLarge_4	30503	10835335

1.2.1.1 Coil: FL1

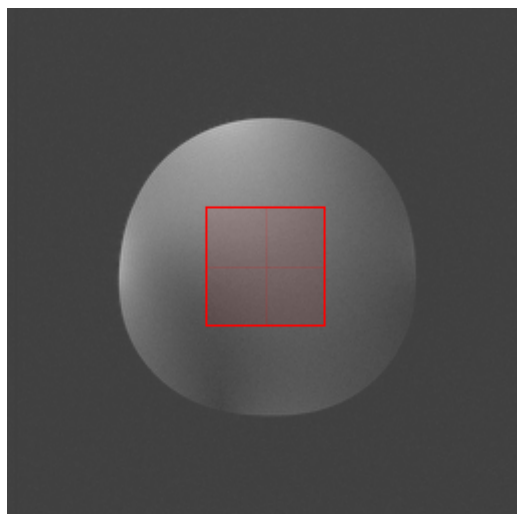

Signal

S/N	
Parabola Fit	Tolerance: ≥ 25.0
28.1	

Brightness	
Parabola Fit	Tolerance: 350 ... 1,400
509	


Homogeneity and Transmitter Amplitude	
Variable	[V]

(continued)

Transmitter Amplitude

385

1.2.1.2 Coil: FL2

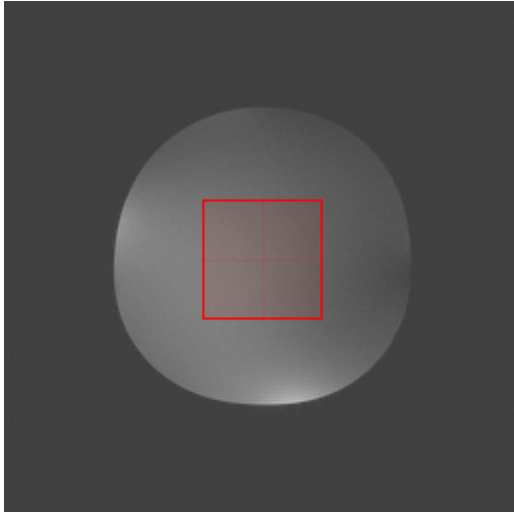

Signal

S/N

Parabola Fit

Tolerance: ≥ 25.0

38.1


Brightness

Parabola Fit

814

Homogeneity and Transmitter Amplitude

Variable	[V]
Transmitter Amplitude	385

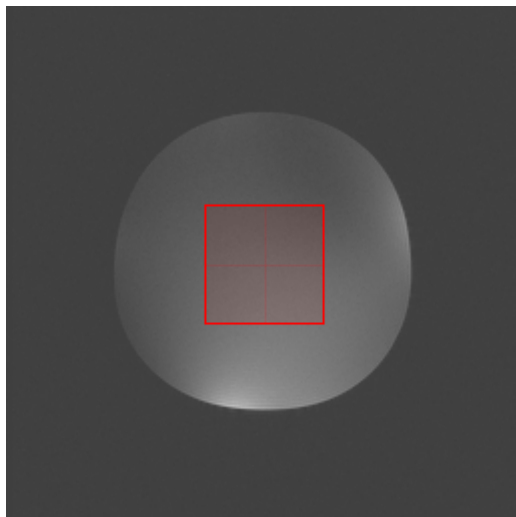1.2.1.3 Coil: FL3

Signal

S/N

Parabola Fit


Tolerance: ≥ 25.0 **34.2****Brightness**

Parabola Fit

665

Homogeneity and Transmitter Amplitude

Variable	[V]
Transmitter Amplitude	385

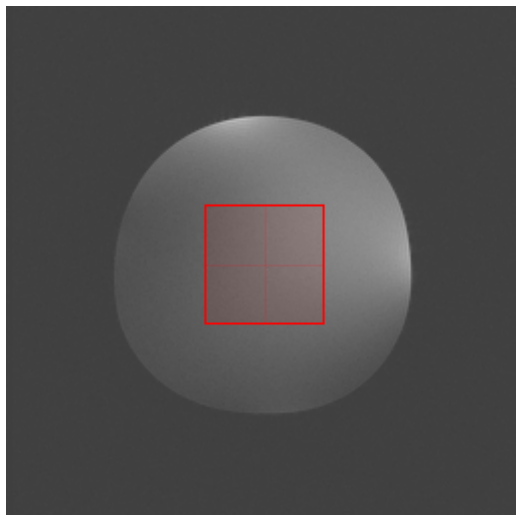1.2.1.4 Coil: FL4

Signal

S/N

Parabola Fit

Tolerance: ≥ 25.0 **37.2****Brightness**

Parabola Fit

749

Homogeneity and Transmitter Amplitude	
Variable	[V]
Transmitter Amplitude	385

Siemens Healthcare Headquarters

Siemens Healthcare GmbH

Henkestr. 127

91052 Erlangen

Germany

Phone: +49 9131 84-0

[siemens.com/healthcare](https://www.siemens.com/healthcare)