


MAGNETOM Vida

Customer Quality Assurance

Serial Number: 175692

Customer Quality Assurance

Demographical Data					
Customer		System		Service	
Customer Name		Product Name	MAGNETOM Vida	Service Center Name	Customer Care Center
Customer ID	400-567654	Material Number	11060815	Street	Wheeler Road
Hospital Name	Beth Israel Deaconess	Serial Number	175692	Street Number	200 suite 303
Street	Brookline Ave	Equipment Number		Zip Code	
Street Number	330	Modality	MR	City	Boston
Zip Code	02215	Handed Over Date	6/28/2018	District	
City	Boston	Software Version	syngo MR XA11	Country	US
District	NorthEast			Phone	18008887436
Country	US			Fax	
Customer Phone Number	16179759913			E-Mail	uptime-services.team@siemens-healthineers.com

Summary		
Function	Start Time	Result
Coil Check - Spine 32 RS (1100)	5/14/2020 2:59:22 PM	 Not OK

1 Coil Check - Spine 32 RS (1100)

1.1 Summary


Summary	
	Status
ExamDb:SIEMENS/System/CoilQA/Spine_32_RS/025_Spine_32_RS_SP1_Tra	OK
ExamDb:SIEMENS/System/CoilQA/Spine_32_RS/025_Spine_32_RS_SP2_Tra	✖ Not OK
ExamDb:SIEMENS/System/CoilQA/Spine_32_RS/025_Spine_32_RS_SP3_Tra	OK
ExamDb:SIEMENS/System/CoilQA/Spine_32_RS/025_Spine_32_RS_SP4_Tra	OK
ExamDb:SIEMENS/System/CoilQA/Spine_32_RS/025_Spine_32_RS_SP5_Tra	OK
ExamDb:SIEMENS/System/CoilQA/Spine_32_RS/025_Spine_32_RS_SP6_Tra	OK
ExamDb:SIEMENS/System/CoilQA/Spine_32_RS/025_Spine_32_RS_SP7_Tra	✖ Not OK
ExamDb:SIEMENS/System/CoilQA/Spine_32_RS/025_Spine_32_RS_SP8_Tra	OK
Program Result	✖ Out of Spec

1.2 Details

1.2.1 Protocol: ExamDb:SIEMENS/System/CoilQA/Spine_32_RS/025_Spine_32_RS_SP1_Tra

Coil Info		
Coil ID	Serial Number	Part Number
Spine_32_RS	1100	10835689

1.2.1.1 Coil: S11

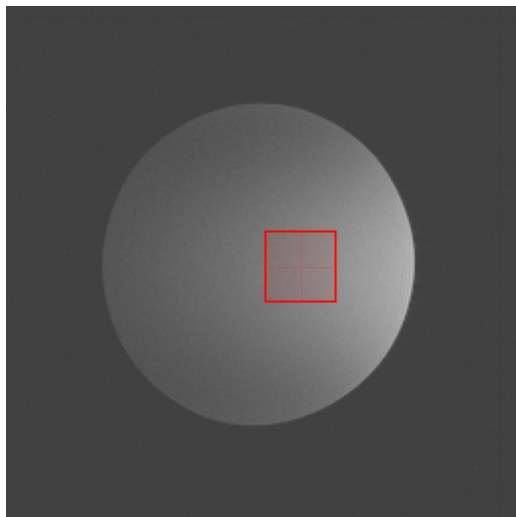


Signal

S/N
Parabola Fit
Tolerance: ≥ 30.0
47.9

Brightness
Parabola Fit
Tolerance: 300 ... 1,200
518

Homogeneity and Transmitter Amplitude

Variable	[V]
Transmitter Amplitude	441

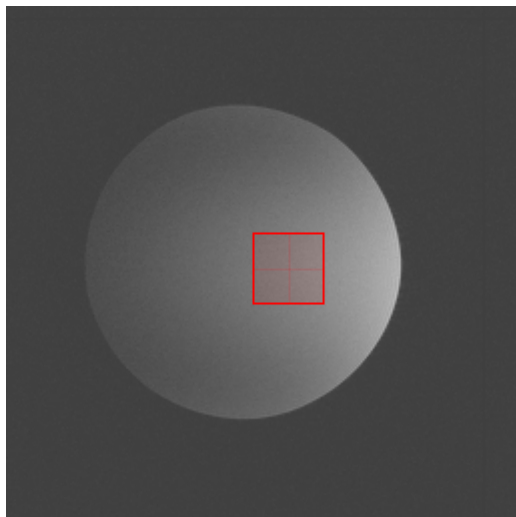1.2.1.2 Coil: S12

Signal

S/N

Parabola Fit


Tolerance: ≥ 30.0 **41.2****Brightness**

Parabola Fit

475

Homogeneity and Transmitter Amplitude

Variable	[V]
Transmitter Amplitude	441

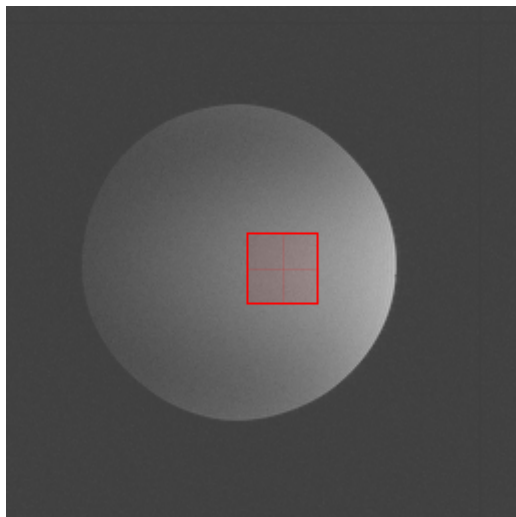1.2.1.3 Coil: S13

Signal

S/N

Parabola Fit


Tolerance: ≥ 30.0 **36.8****Brightness**

Parabola Fit

371

Homogeneity and Transmitter Amplitude	
Variable	[V]
Transmitter Amplitude	441

1.2.1.4 Coil: S14

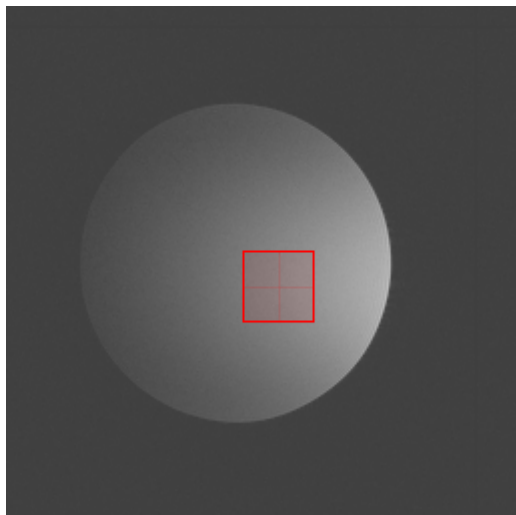

Signal

S/N
Parabola Fit
Tolerance: ≥ 46.0
57.0


Brightness
Parabola Fit
781

Homogeneity and Transmitter Amplitude	
Variable	[V]
Transmitter Amplitude	441

1.2.2 Protocol: ExamDb:SIEMENS/System/CoilQA/Spine_32_RS/025_Spine_32_RS_SP2_Tra

Coil Info		
Coil ID	Serial Number	Part Number
Spine_32_RS	1100	10835689

1.2.2.1 Coil: S21

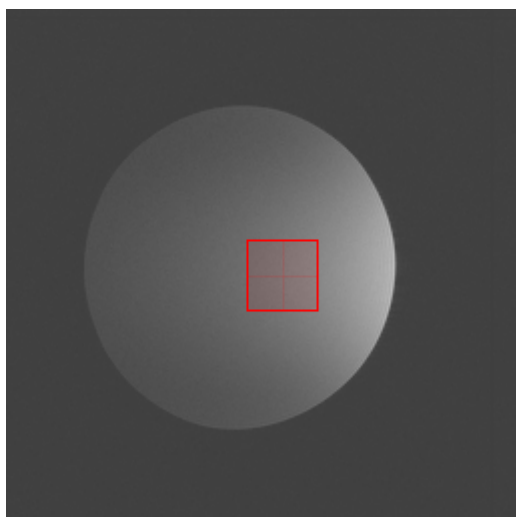


Signal

S/N
Parabola Fit
Tolerance: ≥ 32.0
50.2

Brightness
Parabola Fit
679

Homogeneity and Transmitter Amplitude	
Variable	[V]
Transmitter Amplitude	433

1.2.2.2 Coil: S22

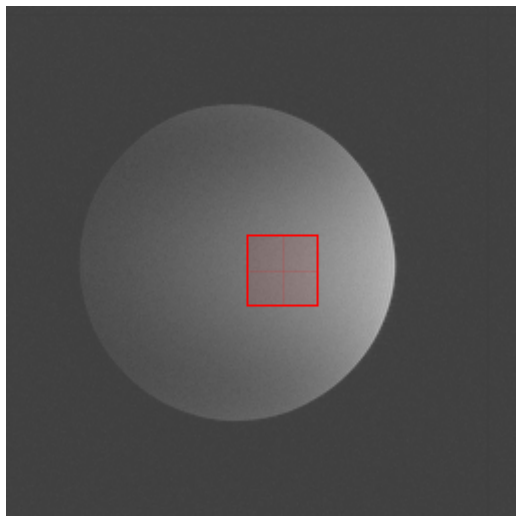

Signal


S/N
Parabola Fit
37.6

Tolerance: ≥ 29.0

Brightness
Parabola Fit
407

Homogeneity and Transmitter Amplitude	
Variable	[V]
Transmitter Amplitude	433

1.2.2.3 Coil: S23

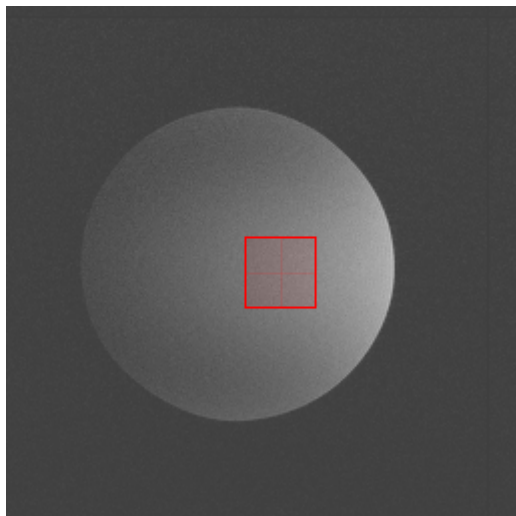

Signal


S/N
Parabola Fit
Tolerance: ≥ 27.0
✖ 23.8



Brightness
Parabola Fit
388

Homogeneity and Transmitter Amplitude	
Variable	[V]
Transmitter Amplitude	433

1.2.2.4 Coil: S24

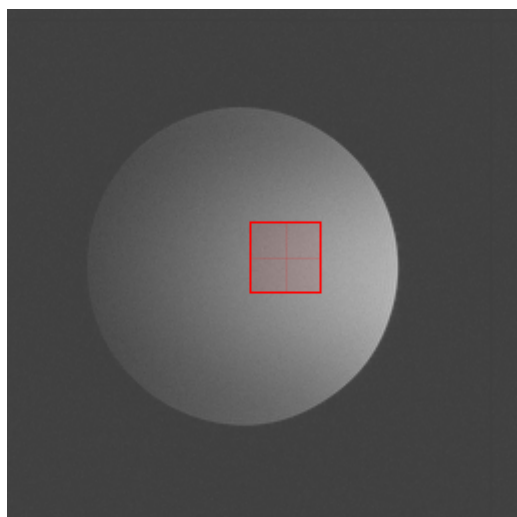

Signal

S/N
Parabola Fit
Tolerance: ≥ 37.0
48.1


Brightness
Parabola Fit
511

Homogeneity and Transmitter Amplitude	
Variable	[V]
Transmitter Amplitude	433

1.2.3 Protocol: ExamDb:SIEMENS/System/CoilQA/Spine_32_RS/025_Spine_32_RS_SP3_Tra

Coil Info		
Coil ID	Serial Number	Part Number
Spine_32_RS	1100	10835689

1.2.3.1 Coil: S31

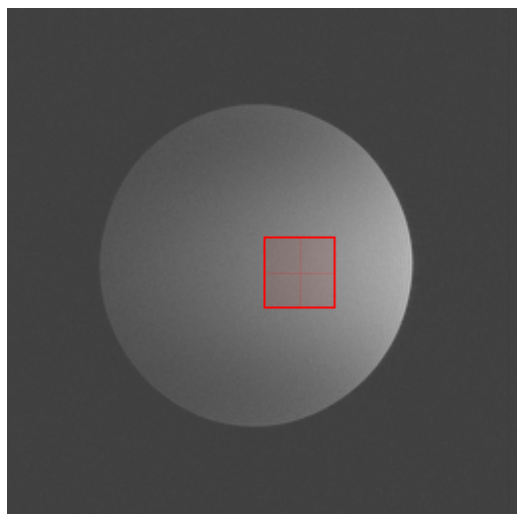


Signal

S/N	
Parabola Fit	Tolerance: ≥ 30.0
48.4	

Brightness	
Parabola Fit	
557	

Homogeneity and Transmitter Amplitude	
Variable	[V]
Transmitter Amplitude	422

1.2.3.2 Coil: S32

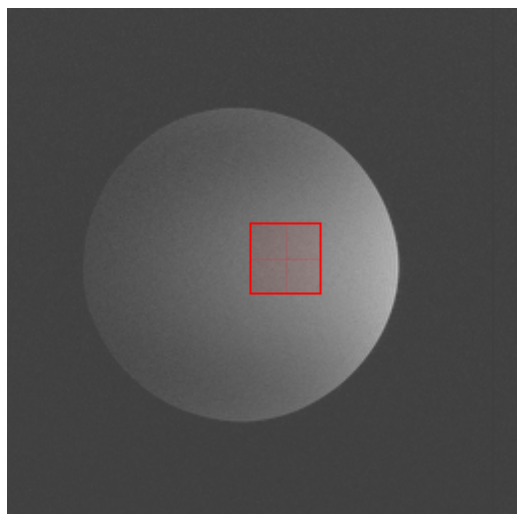


Signal

S/N	
Parabola Fit	Tolerance: ≥ 30.0
30.5	

Brightness	
Parabola Fit	
519	

Homogeneity and Transmitter Amplitude	
Variable	[V]
Transmitter Amplitude	422

1.2.3.3 Coil: S33

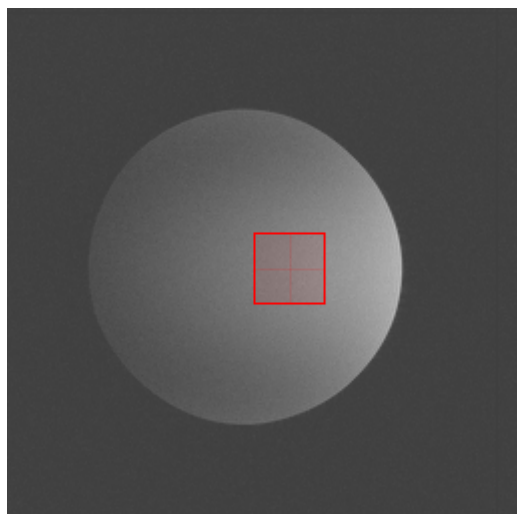


Signal

S/N	
Parabola Fit	Tolerance: ≥ 30.0
37.6	

Brightness	
Parabola Fit	
394	

Homogeneity and Transmitter Amplitude	
Variable	[V]
Transmitter Amplitude	422

1.2.3.4 Coil: S34

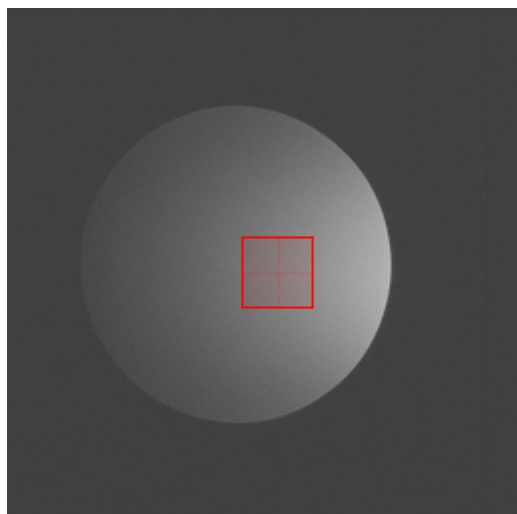

Signal

S/N	
Parabola Fit	Tolerance: ≥ 46.0
49.2	


Brightness	
Parabola Fit	
710	

Homogeneity and Transmitter Amplitude	
Variable	[V]
Transmitter Amplitude	422

1.2.4 Protocol:
ExamDb:SIEMENS/System/CoilQA/Spine_32_RS/025_Spine_32_RS_SP4_Tra

Coil Info		
Coil ID	Serial Number	Part Number
Spine_32_RS	1100	10835689

1.2.4.1 Coil: S41

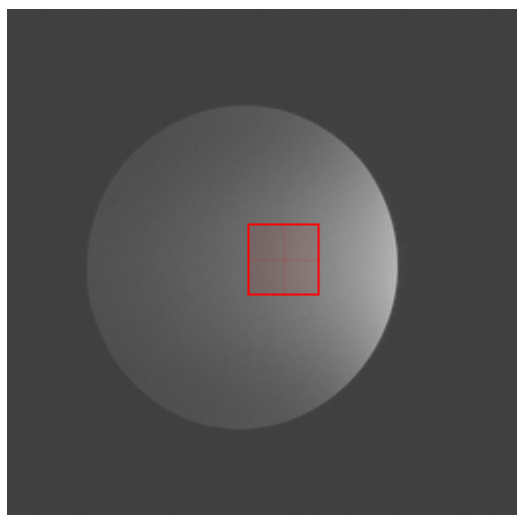

Signal

S/N
Parabola Fit
Tolerance: ≥ 32.0
51.2

Brightness


Parabola Fit

682**Homogeneity and Transmitter Amplitude**

Variable

[V]

Transmitter Amplitude

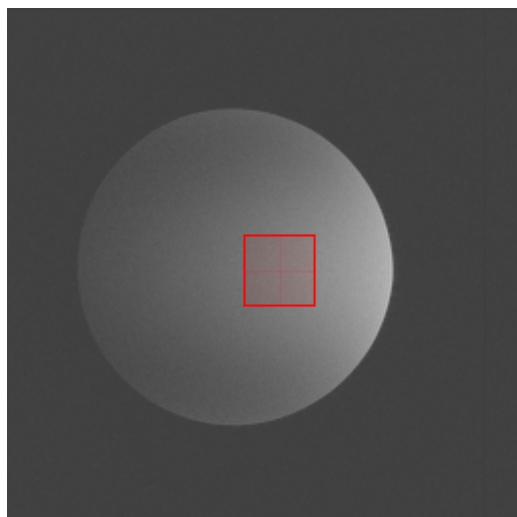422**1.2.4.2 Coil: S42**

Signal

S/N

Parabola Fit

Tolerance: ≥ 29.0 **38.0**

Brightness


Parabola Fit

394**Homogeneity and Transmitter Amplitude**

Variable

[V]

Transmitter Amplitude

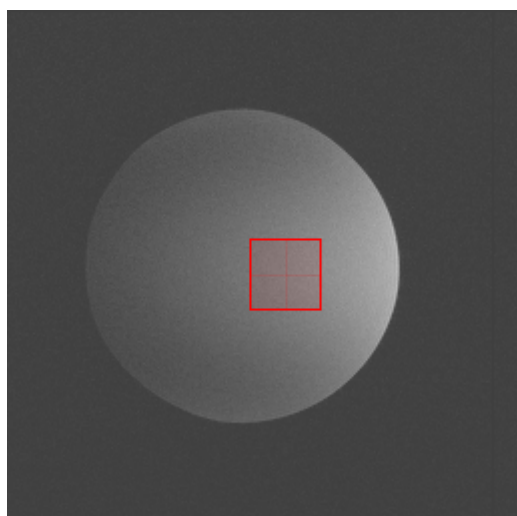422**1.2.4.3 Coil: S43**

Signal

S/N

Parabola Fit

Tolerance: ≥ 27.0 **28.4**

Brightness


Parabola Fit

393**Homogeneity and Transmitter Amplitude**

Variable

[V]

Transmitter Amplitude

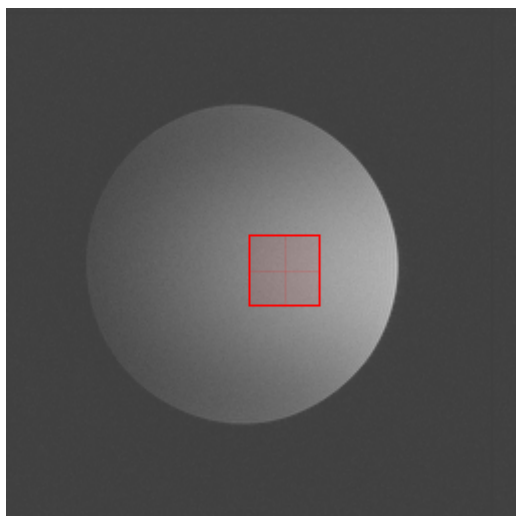422**1.2.4.4 Coil: S44**

Signal

S/N

Parabola Fit

Tolerance: ≥ 37.0 **47.1**


Brightness
Parabola Fit
549

Homogeneity and Transmitter Amplitude	
Variable	[V]
Transmitter Amplitude	422

1.2.5 Protocol:
ExamDb:SIEMENS/System/CoilQA/Spine_32_RS/025_Spine_32_RS_SP5_Tra

Coil Info		
Coil ID	Serial Number	Part Number
Spine_32_RS	1100	10835689

1.2.5.1 Coil: S51

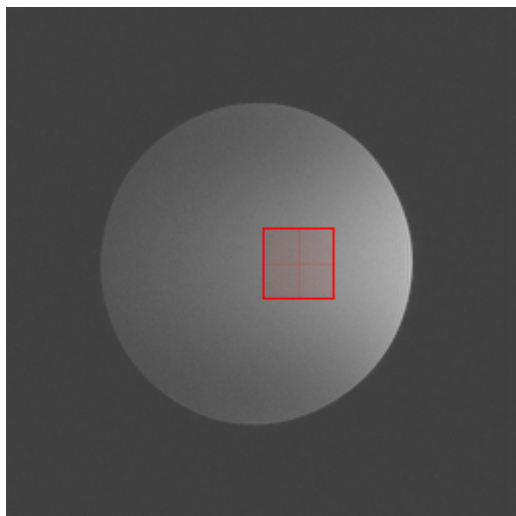


Signal

S/N
Parabola Fit
Tolerance: ≥ 30.0
48.0

Brightness
Parabola Fit
523

Homogeneity and Transmitter Amplitude	
Variable	[V]
Transmitter Amplitude	420

1.2.5.2 Coil: S52

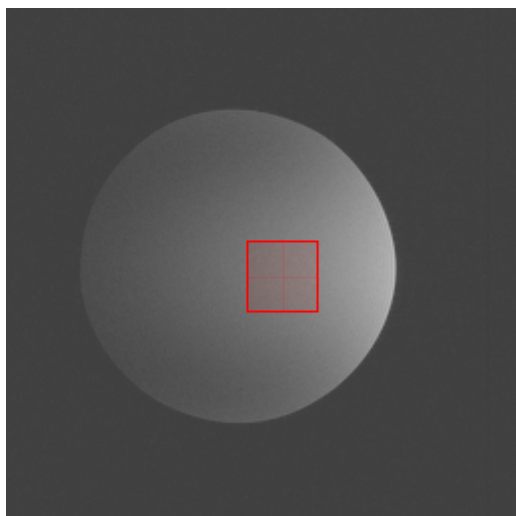


Signal

S/N
Parabola Fit
Tolerance: ≥ 30.0
42.0

Brightness
Parabola Fit
473

Homogeneity and Transmitter Amplitude	
Variable	[V]
Transmitter Amplitude	420

1.2.5.3 Coil: S53

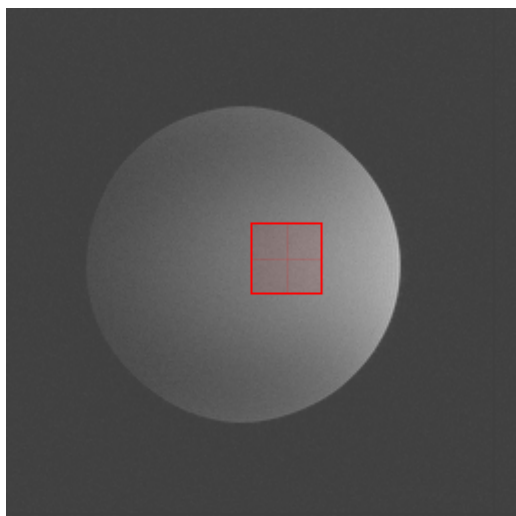


Signal

S/N
Parabola Fit
Tolerance: ≥ 30.0
34.8

Brightness
Parabola Fit
390

Homogeneity and Transmitter Amplitude	
Variable	[V]
Transmitter Amplitude	420

1.2.5.4 Coil: S54

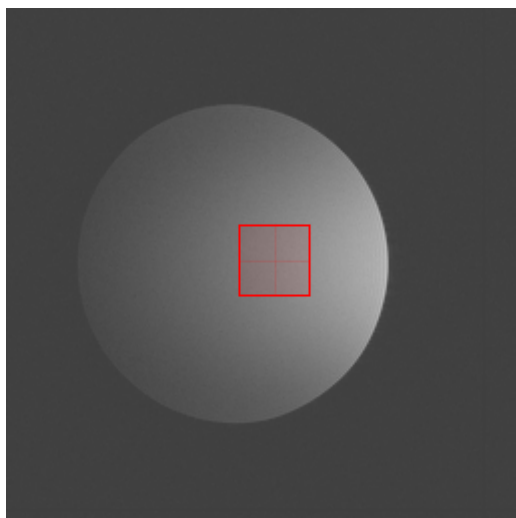

Signal

S/N
Parabola Fit
Tolerance: ≥ 46.0
58.2


Brightness
Parabola Fit
734

Homogeneity and Transmitter Amplitude	
Variable	[V]
Transmitter Amplitude	420

1.2.6 Protocol:
ExamDb:SIEMENS/System/CoilQA/Spine_32_RS/025_Spine_32_RS_SP6_Tra

Coil Info		
Coil ID	Serial Number	Part Number
Spine_32_RS	1100	10835689

1.2.6.1 Coil: S61

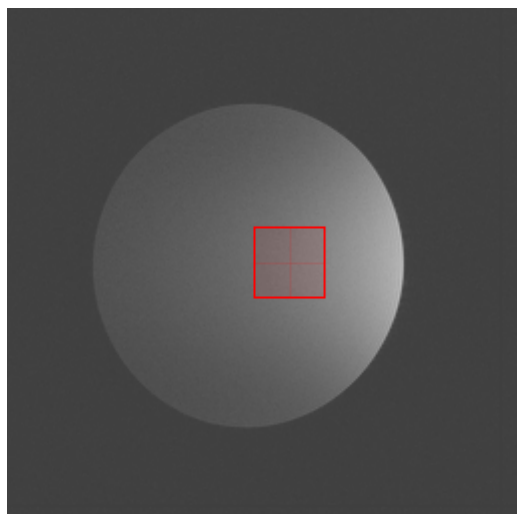


Signal

S/N	
Parabola Fit	Tolerance: ≥ 32.0
43.8	

Brightness	
Parabola Fit	
619	

Homogeneity and Transmitter Amplitude	
Variable	[V]
Transmitter Amplitude	421

1.2.6.2 Coil: S62

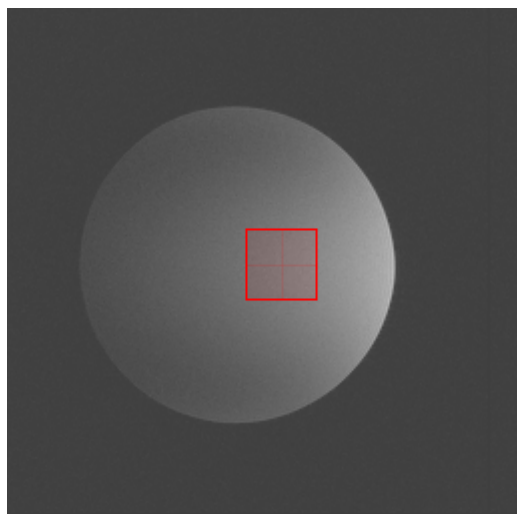


Signal

S/N	
Parabola Fit	Tolerance: ≥ 29.0
38.5	

Brightness	
Parabola Fit	
399	

Homogeneity and Transmitter Amplitude	
Variable	[V]
Transmitter Amplitude	421

1.2.6.3 Coil: S63

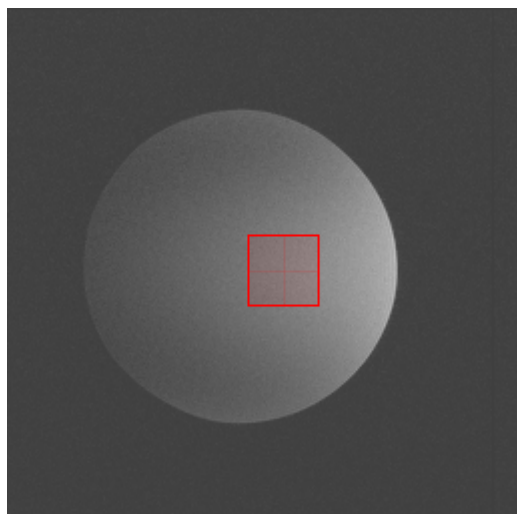


Signal

S/N	
Parabola Fit	Tolerance: ≥ 27.0
28.0	

Brightness	
Parabola Fit	
377	

Homogeneity and Transmitter Amplitude	
Variable	[V]
Transmitter Amplitude	421

1.2.6.4 Coil: S64

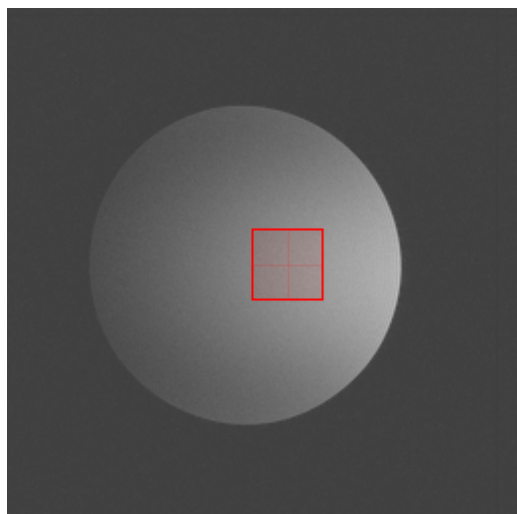

Signal

S/N	
Parabola Fit	Tolerance: ≥ 37.0
47.0	


Brightness	
Parabola Fit	
514	

Homogeneity and Transmitter Amplitude	
Variable	[V]
Transmitter Amplitude	421

1.2.7 Protocol:
ExamDb:SIEMENS/System/CoilQA/Spine_32_RS/025_Spine_32_RS_SP7_Tra

Coil Info		
Coil ID	Serial Number	Part Number
Spine_32_RS	1100	10835689

1.2.7.1 Coil: S71

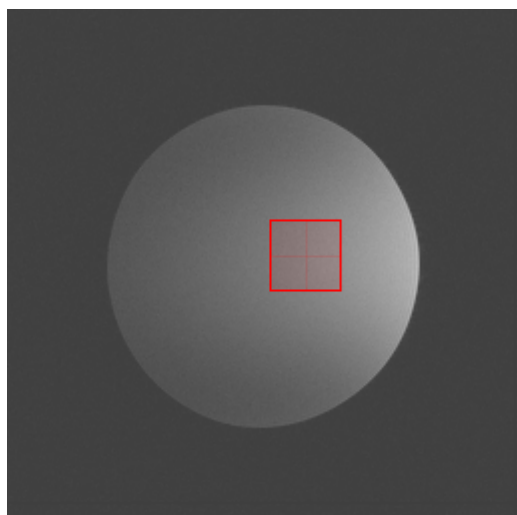

Signal

S/N
Parabola Fit
Tolerance: ≥ 30.0
45.7

Brightness


Parabola Fit

523**Homogeneity and Transmitter Amplitude**

Variable

[V]

Transmitter Amplitude

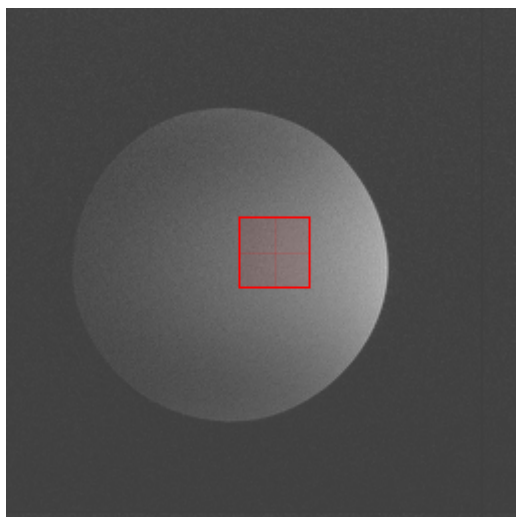434**1.2.7.2 Coil: S72**

Signal

S/N


Parabola Fit

Tolerance: ≥ 30.0 **✖ 22.9**

Brightness
Parabola Fit
430

Homogeneity and Transmitter Amplitude	
Variable	[V]
Transmitter Amplitude	434

1.2.7.3 Coil: S73

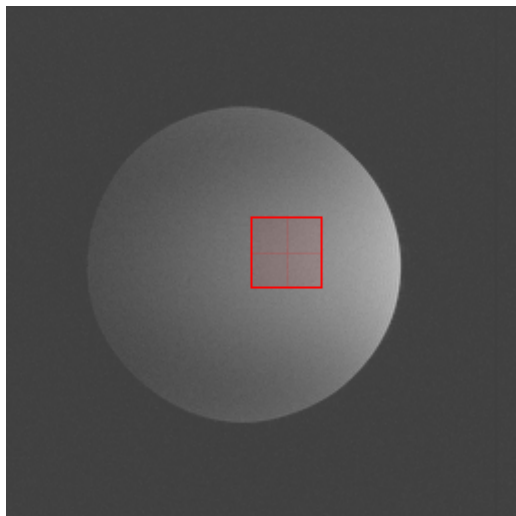


Signal

S/N
Parabola Fit
Tolerance: ≥ 30.0
37.5

Brightness
Parabola Fit
386

Homogeneity and Transmitter Amplitude	
Variable	[V]
Transmitter Amplitude	434

1.2.7.4 Coil: S74

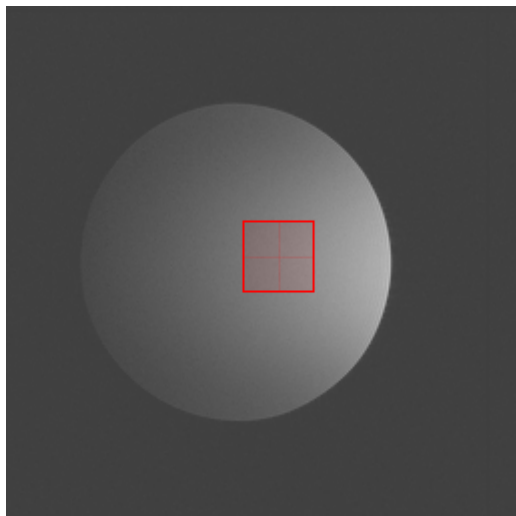

Signal

S/N
Parabola Fit
Tolerance: ≥ 46.0
55.5

Brightness
Parabola Fit
757

Homogeneity and Transmitter Amplitude	
Variable	[V]
Transmitter Amplitude	434

1.2.8 Protocol: ExamDb:SIEMENS/System/CoilQA/Spine_32_RS/025_Spine_32_RS_SP8_Tra


Respiratory Sensor Level Check				
Sensor	Type	Condition	Result	Tolerance
0	Signal	Loaded	4,393,952	4,000,000 ... 5,500,000
0	Signal	Unloaded	4,415,402	4,000,000 ... 5,500,000
0	Signal Difference	Unloaded - Loaded	21,450	10,000 ... 250,000
0	Noise	Loaded	65	≤ 100
0	Noise	Unloaded	68	≤ 100
1	Signal	Loaded	4,639,487	4,000,000 ... 5,500,000
1	Signal	Unloaded	4,696,209	4,000,000 ... 5,500,000

(continued)

Respiratory Sensor Level Check				
Sensor	Type	Condition	Result	Tolerance
1	Signal Difference	Unloaded - Loaded	56,723	10,000 ... 250,000
1	Noise	Loaded	55	≤ 100
1	Noise	Unloaded	75	≤ 100

Coil Info		
Coil ID	Serial Number	Part Number
Spine_32_RS	1100	10835689

1.2.8.1 Coil: S81

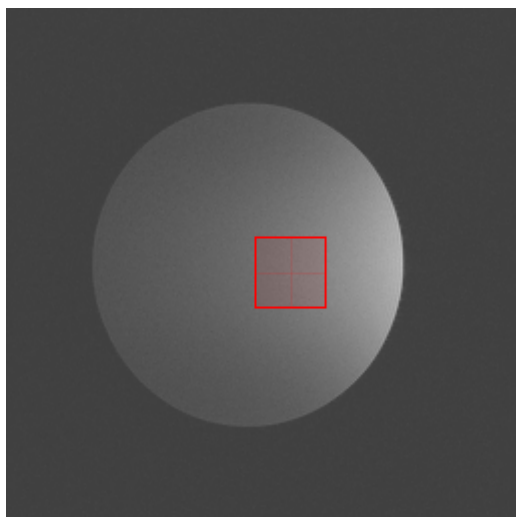


Signal

S/N
Parabola Fit
Tolerance: ≥ 35.0
36.5

Brightness
Parabola Fit
599

Homogeneity and Transmitter Amplitude	
Variable	[V]
Transmitter Amplitude	455

1.2.8.2 Coil: S82

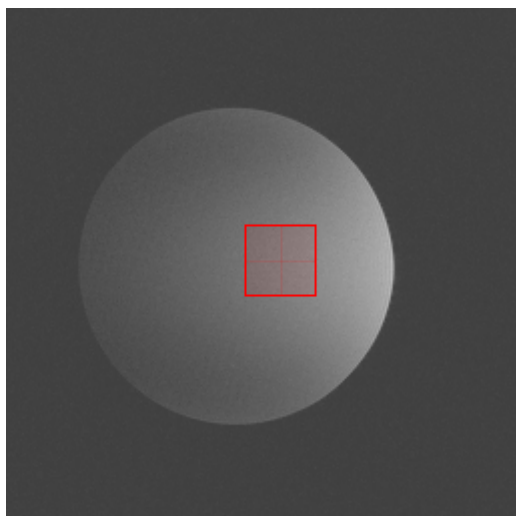


Signal

S/N
Parabola Fit
Tolerance: ≥ 25.0
32.8

Brightness
Parabola Fit
388

Homogeneity and Transmitter Amplitude	
Variable	[V]
Transmitter Amplitude	455

1.2.8.3 Coil: S83

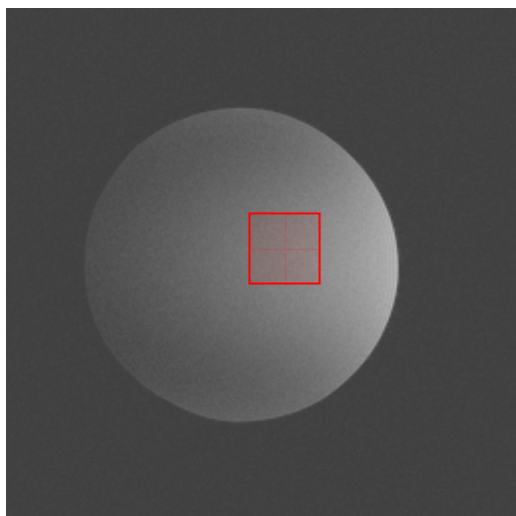


Signal

S/N
Parabola Fit
Tolerance: ≥ 25.0
28.3

Brightness
Parabola Fit
359

Homogeneity and Transmitter Amplitude	
Variable	[V]
Transmitter Amplitude	455

1.2.8.4 Coil: S84

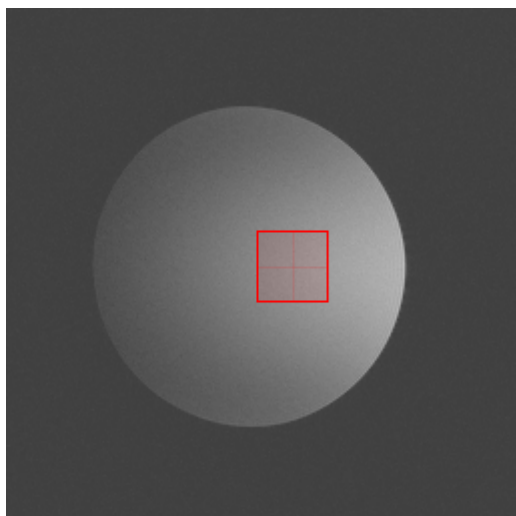

Signal

S/N	
Parabola Fit	Tolerance: ≥ 35.0
47.0	

Brightness	
Parabola Fit	
535	

Homogeneity and Transmitter Amplitude	
Variable	[V]
Transmitter Amplitude	455

Siemens Healthcare Headquarters

Siemens Healthcare GmbH

Henkestr. 127

91052 Erlangen

Germany

Phone: +49 9131 84-0

[siemens.com/healthcare](https://www.siemens.com/healthcare)